

## What does good value mean to you?

"I want value as I define it: quality equipment that helps me care for my patients. Ritter exam tables provide exactly what I need. They're not only cost effective, but provide years of durability. Ritter exam tables make it easier for me and my staff."

### Ritter 204

## S P E C I F I C A T I O N S



*Seamless UltraFree Upholstery Top*  
The Seamless UltraFree upholstery top provides a comfortable and easy-to-clean option for your practice, with PVC-free material.

### Ritter 204 Manual Examination Table

#### Specifications

##### Height with Upholstery:

Soft Touch Top -  
32.13" (81.6 cm)  
Seamless Top -  
33.13" (84.2 cm)

Patient Weight Capacity: 500 lb  
(226.8 kg)

##### Electrical Receptacle

(115 VAC 5 A Max.)  
(230 VAC 2.5 A Max.)  
Duplex, hospital grade  
electrical receptacle on  
left side

##### Soft Touch Upholstery Top:

30.5" W x 56" L  
(77.5 x 142.5 cm)

Paper Roll Holder -  
Holds 21" x 4.75" Paper  
Roll (53.3 cm x 12 cm)

Length with footrest  
extended: 73.75"  
(187.3 cm)

##### Seamless Upholstery Top:

28.5" W x 59.5" L  
(72.4 x 151.1 cm)

Paper Roll Holder -  
Holds 21" x 3.5" Paper  
Roll (53.3 cm x 8.9 cm)

Length with Footrest  
extended: 76" (193 cm)

UltraFree Upholstery Top:  
Specs to come...

#### Options

Soft Touch, Seamless  
or Seamless UltraFree  
Upholstery Options

Pelvic tilt  
(1.5 degrees or 7 degrees)

Drawer heater (45 W)\*

#### Accessories

Knee crutches

Lower body drain pan

Stainless steel pan

Assist handles

Adjustable arm system

Flat orthopedic upholstery

#### Electrical

115 VAC, 50/60 Hz, 5.0 A  
(with receptacle)

115 VAC, 50/60 Hz, 5.5 A  
(with drawer heater)

230 VAC, 50/60 Hz, 2.5 A  
(with receptacle)

UL 60601-1 and Can/CSA  
601.1-M90 Standard for  
Safety

\* Not available for  
sale in Europe.

#### Standard Upholstery Colors

- Mist-813
- Stone-814
- Restful Path-858
- Iron Ore-851
- Lunar Gray-845
- Obsidian-857
- Oasis-843
- Soothing Blue-855
- Dream-847
- Mineral-844
- Healing Waters-848
- Shaded Garden-853
- Citrus-849
- Sandy Retreat-850
- Curative Copper-846
- Cranberry-859
- Dark Linen-856
- Latte-860
- Robust Brown-852
- Deep Earth-854

#### UltraFree\* Upholstery Colors

- UltraFree Stone-865
- UltraFree Restful Path-868
- UltraFree Obsidian-867
- UltraFree Harbor-861
- UltraFree Wheatgrass-862
- UltraFree Spice-863
- UltraFree Cranberry-869
- UltraFree Dark Linen-866
- UltraFree Latte-870
- UltraFree Branch-864

#### \*Ultraleather® by Ultrafabrics

## CE

### CARB 93120.2 Phase 2 Compliant

Midmark is an ISO 13485 and  
ISO 9001 Certified Company.

For more information  
or a demonstration, contact  
your Midmark dealer or Midmark  
Corporation, Dayton, OH.  
Call: 1-800-MIDMARK  
Fax: 1-800-365-8631  
midmark.com

Authorized CE Representative  
Countries in the EEC should direct all questions, incidents and  
complaints to Midmark's Authorized CE representative listed below:

Midmark EMEA Ltd  
Beech House, First Floor  
East Wing, Ancells Business Park  
Fleet, Hampshire,  
GU51 2UN, UK  
Call: +44 (0) 1252 360 940  
Fax: +44 (0) 1252 360 941

© Midmark Corporation 2005  
Midmark Corporation, Dayton, OH.  
Products subject to improvement changes without notice.  
Litho in U.S.A. 007-0560-00 Rev. J1 (1/18)

**Ritter**  
by MIDMARK



**Ritter**  
by MIDMARK

# 204 Manual Examination Table



## We asked the experts

What makes Ritter the most purchased brand of exam tables in the industry? It's our willingness to research, study and listen so we can meet your needs—and the needs of your patients.

You need...

- three side drawers for optimal supply management
- modern styling to coordinate with your office
- an updated, ergonomic top design to help patients feel comfortable
- additional weight capacity up to 500 pounds

So, what was our solution? The Ritter 204 Manual Exam Table—the latest in design, comfort and functionality in a manual exam table.



**One-Piece Seamless Drawers**  
Designed with easy-to-clean, rounded edges that are seamless to contain spilled fluids and help improve infection control.


**Spacious Side Drawers** - Three spacious side drawers provide plenty of storage volume for gowns, drapes and other items, yet are made to help organize your assortment of smaller supplies and bulk supplies. In addition, there are two front drawers for easy access to items during pelvic exams.



**Stylish and Easy to Clean** - The seamless design of the upholstered top is as soft and comfortable as it looks, leaving your patients feeling more relaxed. And the seamless upholstery is designed to eliminate areas where dirt and bacteria can accumulate so clean-up is easier, faster and more effective.



**Removable Top System** - The one-piece, upholstered top is designed to be easily removed without tools, allowing for thorough cleaning or a quick upholstery color change.

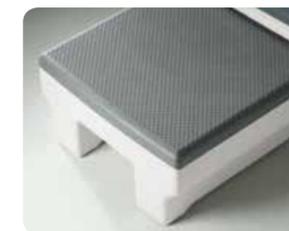

**Enlarged Footstep** - We increased the size of the footstep over 40 percent above the industry standard. It is designed to allow patients to feel more secure and get on the table more confidently.



**Storage Access and Volume** - Our pass-through side drawers provide access to supplies from either side of the table and give you the ability to change your room configuration.


## Accessories that add convenience



**Assist Handles**  
Designed to provide patients a secure means of access to the exam table or to reposition themselves during an exam. Easily removed without tools, assist handles offer greater patient access.




**Flat Upholstery**  
Provides an excellent working surface for orthopedic, chiropractic or osteopathic uses.

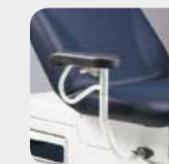


**Adjustable Arm System**  
Designed to offer patients stable support for accessing the exam table and a sense of security during transfer. Arms rotate and lock into one of three positions and are easily removable.



**Ritter 253 LED Exam Light**  
Designed to be the perfect complement to any procedure room, the 253 provides features like easy access controls, high intensity light and an adjustable focal spot. Available in three mounting options.

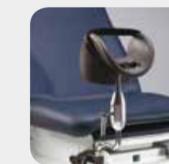

**Articulating Knee Crutches**  
A more ergonomic and comfortable option compared to standard stirrups. The ball joint design offers ultimate flexibility in patient positioning.

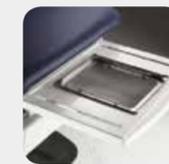

**Lower Body Drain Pan**  
Ideal for lower body procedures or other situations in which removal of fluid and debris is necessary.



**IQvitals® Table Mount**  
Convenient table mount for the Midmark IQvitals vital signs device, designed to fit all Midmark and Ritter exam tables.



**Standard Knee Crutches**  
Provide leg support for extended procedures where fatigue is a factor or when the patient lacks sufficient leg strength or control. Swivel to adjust for knee and leg support.



## Everyone is entitled to equal access care.

To really care for people, you must understand the patients you care for—today and tomorrow. Due to an increasingly aging population, your patients may need assistance accessing the exam table. With a low 18" height, Barrier-Free® Exam Tables are designed to make it easy for you to accommodate all patients with dignity, regardless of mobility concerns. Help them help themselves.